


Company Introduction

Beijing Globalipl Development CO., Ltd. is a high-tech company engaged in development, research, producing, sale and after-sale of beauty equipment. Nowadays, we provide OEM service for several big national suppliers on fittings, such as PCB and probes. Our main business focuses on developing and production of high quality beauty equipment, which create our leadership in this field.



With high technology, stable quality and perfect treatment result, our equipment has gained good reputation from numerous customers. Furthermore, our company has achieved great competitiveness in the world market, simultaneously, promoted our technology development in a high speed.

Our company insists on people oriented, front-running sciences & tech and client first principle, taking survival by high quality and development through sciences & tech. We provide perfect after-sale service to make customer satisfied, and high performance-price ratio products to make customer get more profit. OEM Service is our speciality. We can produce software and hardware according to your requirements. Make it easy for you to become a professional manufacturer!

Sincerely wish we have a chance to cooperate with each other in the near future!

Our mission : To serve people who love beauty around the world.

Our values: 1.Individual and team: we enjoy the honour and accept the failure together 2.Technique is the base of the company 3.Quality is the life of the company 4.Focusing on our major product 5.Creating world class brand-Globalipl.

Product Catalogue

01 SHR IPL Machine



US601F US001F US610 US608F

16 IPL / SHR + Diode Laser Machine



US419

06 DPL-Dye Pulse Light Machine



US002F

18 4H (IPL+RF+E-Light+YAG Laser) Machine



US002H US609H

08 Diode Laser Machine



US450 US415 US417

20 3H (IPL+RF+E-Light) Machine



US606S

12 Micro Channel Diode Laser Machine



US460 US420

22 IPL E-Light Machine



US002 US609B

14 Fiber-Coupled Diode Laser Machine



US408G US417G

24 Cryotherapy Slimming Machine



US09S US320C US601N

Product Catalogue

28

Slimming Machine—RF Vacuum Cavitation Diode laser



US07B



US307



US308

40

Q Switch ND YAG Laser Machine



US409



US406



US400N

30

Vacuum Roller Slimming Machine



US08S

45

Laser Vascular Removal Machine



US410C

32

CO₂ Fractional Laser Machine



US800



US900



US800N



US900N

46

H₂O₂ Facial Jet Peel Machine



US710A



US710HM

36

PicoSecond Laser Machine



US506



US508

48

Air Cooling Machine



Cool -01

38

Q Switch ND YAG Laser Machine



US412



US412C



US407



US407A

3^{plus} SHR IPL Machine

Flexible Screen, 3 Handpieces

Functions

- Super hair removal
- Super skin rejuvenation
- Pigment therapy
- Wrinkle removal
- Vascular therapy
- Breast liftup



The screen can be folded up and down.



Rotate 180° left and right



US601F

Three handpieces for different parts of body treatment, more convenient



8x15mm² for armpit, and bikini hair removal



12x30mm² for facial treatment



15x50mm² for large area hair removal


Work theory

SHR works by gradually heating the dermis to a temperature that effectively damages the hair follicles and prevents regrowth, while avoiding injury to the surrounding tissue. A high repetition rate of short pulses are delivered deep into the dermis, achieving high average power and therapeutically effective heat build-up, with virtually no pain.



OPT perfect pulse - advantage

- Uniform output energy
- Less pain
- Better treatment results
- The speed of light Shot speed: 10 shots per second, faster

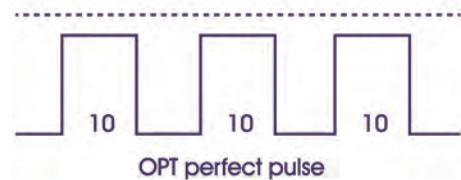

Acne removal




Vascular removal



Pigment removal



Hair removal

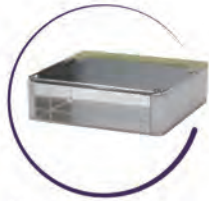

Sapphire

more than 95% light transmittance



Connector

"plug & play", USA & Germany standard



Water pump

no noise, large water flow for better cooling




Xenon lamp

Germany Heraeus, 1,000,000 shots



Filters

480nm,530nm,590nm,640nm,690nm

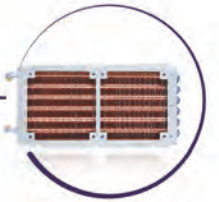

USA Radiator

Imported 316L stainless steel for excellent cooling

**Semiconductor Thermal Management System,
decrease 11 degree within 8 minutes!!! — For option**



Technical parameters

IPL energy	1-50 J/cm ² , adjustable		
Wavelength	480nm,530nm,590nm,640nm,690nm		
Shot frequency	0.1s / 0.2s / 0.3s / 0.5s / 1s / 2s / 3s / 4s		
Pulse mode	Single Pulse		
Cooling system	USA cooler+air+water+cooling gel(Semiconductor Thermal Management system for option)		
Language	English, Italian, Polish, Russian, Turkish, Spanish, French, Portuguese 8 languages in total		
Voltage	220V/110V		
Package	Aluminum alloy case		
Machine Specifications			
Model No.	US601F	US610	US608F
Display	10.4" color touch screen	8.4"inch color touch screen	8.4"color touch screen
Output power	4000w	2200W	2200W
Spot size	8x15mm ² 12x30mm ² 15x50 mm ²	15x50mm ² 12x30mm ²	15x50mm ² 12x30mm ²
Handles number	3 pcs	2 pcs	2 pcs
Machine size	50cm×58cm×110cm	50x50x48cm	41x46x49cm
Gross weight	77kg	35kg	30kg

SHR IPL Machine

Fast & Comfortable & Permanent



US001F

- Super hair removal
- Super skin rejuvenation
- Pigment therapy
- Wrinkle removal
- Vascular therapy
- Breast liftup

SHR IPL Machine

Fast & Comfortable & Permanent



US610

Functions

- Super Hair Removal
- Super skin rejuvenation
- Pigment therapy
- Wrinkle removal
- Vascular therapy
- Breast liftup



US608F

Dye Pulse Light Machine


DPL Technology



- Hair removal: All kinds of hair color, all kinds of skin color.
- Skin rejuvenation: shrink pores, wrinkles removal and whitening, brighten the skin
- Vascular: Congenital Vascular, Vascular after skin resurfacing, red face(Facial flush), Acne red print, Acne rosacea etc.
- Freckle remove: Freckles, Chloasma, Sunburn, Senile plaques, Pigmentation, etc.

Theory:

DPL (Dye Pulse Light) hair removal technology according to the principle of selective pyrolysis of light that a particular wavelength of intense pulsed light irradiation in hair, absorbed by the melanin hair shaft and follicle selective and instantaneous formation of local high temperature, make the hair follicle coagulation necrosis, and does not damage normal skin and sweat glands, inhibit the growth of the hair, the effect of permanent hair removal.

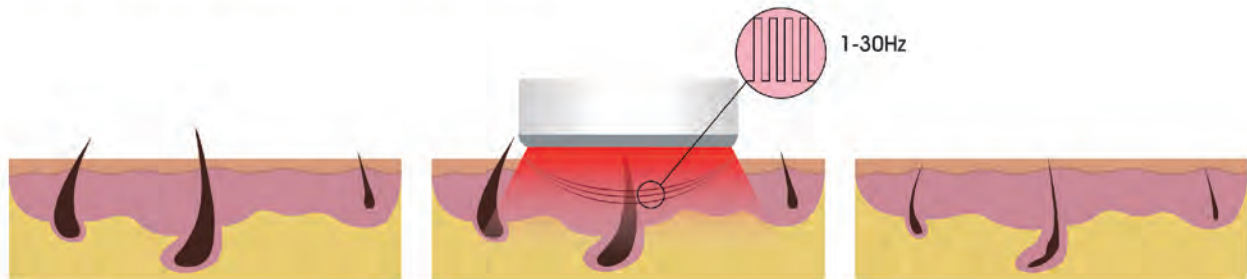

Advantage:

1. Technology upgrade, narrower pulse width and higher frequency.
2. The lesion is more sensitive to DPL light, with higher treatment efficiency and shorter treatment time.
3. Less pain, and improve customer experience.

Handpieces:



- Technology upgrade, narrower pulse width and higher frequency.
- The lesion is more sensitive to DPL light, with higher treatment efficiency and shorter treatment time.
- Less pain, and improve customer experience.

Technical parameters	
Model No.	US002F
Display	10.1" color touch screen
Output Power	4000W
Machine Size	59*52*120cm
Net Weight	60kg
Package	Aluminum Alloy Case
Energy	1-50J/cm ² , adjustable
Frequency	1-30Hz adjustable(DPL); 1-10Hz adjustable(SHR)
Handpiece No.	2 Pcs
Spot Size	15*50mm ² -(DPL) 12*30mm ² -(SHR)
Wavelength	640-950nm-(DPL) 480nm, 530nm, 590nm, 640nm, 690nm-950nm(SHR)
Cooling System	USA cooler+air+water+cooling gel
Skin Cooling	≤0~10 C
Operation Language	English, German, Russian, Turkish, Spanish...16 languages
Voltage	220V/110V

New Sliding Screen Diode Laser Machine

three wavelengths 755nm, 808nm, 1064nm for optional

Function

- Permanent & Painless Hair removal on all body of skin type I-VI.
- lip hair removal, beard hair removal, chest hair removal, armpit hair removal, back hair removal & hair removal on outside bikini line, etc.



US450



US415

Three handpieces for optional



13×39mm²



13×13mm²





13×13mm²

Working Theory

The 755nm wavelength offers more powerful energy absorption by the melanin chromophore, making it ideal for the widest range of hair types and color- especially light-colored and thin hair.

The 808nm wavelength, offers deep penetration of the hair follicle with high average power, a high repetition rate and a large spot size for fast treatment.

The 1064 wavelength is characterized by lower melanin absorption, making it a focused solution for darker skin types.

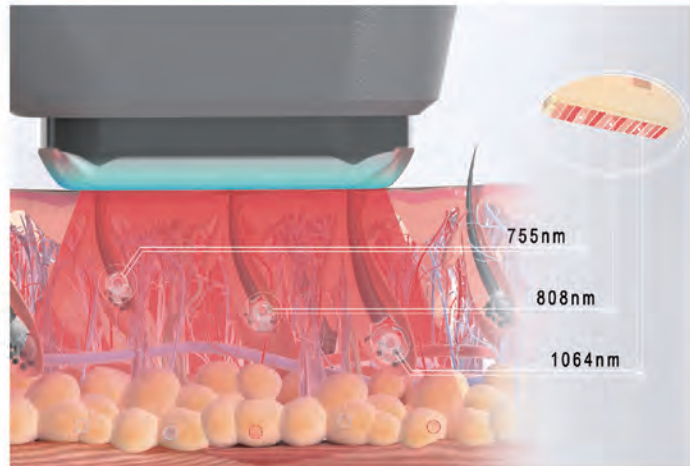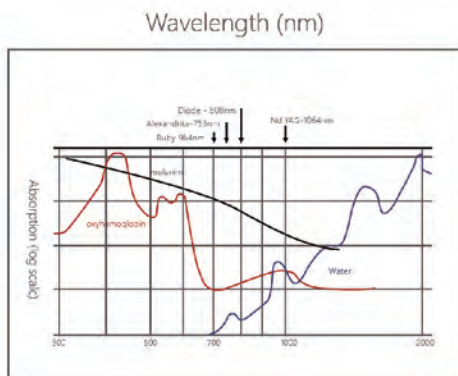

The three wavelengths bars (755nm, 808nm, 1064nm) assembled inside one handle, a light capable of emitting three wavelengths of light, to achieve different skin depths So achieve more thorough hair removal



Light absorbed by melanin



Black pigment disappeared



Hair loss

Palladium bars pure imported from Germany.

- AuSn hard solder assembly technology, higher reliability
- Optimized OE conversion efficiency
- Output power range : 500w-3200w




Germany imported ultra-quiet water pump

- no noise, provide a quiet and comfortable environment for patients
- larger water flow for better cooling, prolong laser and machine life.

New designed integrated handpiece

- Most advanced plug and play technology, no leakage of water and electricity, safe and reliable.
- USA import water connector and Germany import electronic connector.
- Humanized bending design, 100% import tube, maximum protect inner cable from vertical fold, and water flow smoothly for better cooling.

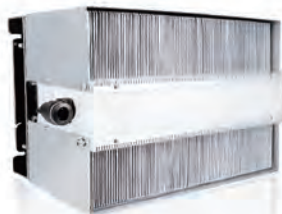

Most advanced semiconductor thermal management cooling system

- The machine could work in a constant temperature environment around 25-28°C, prolong the laser and machine life time, save cost.
- Japanese Semiconductor Thermal Management System ,no noise.



drop 11 degrees
within 8 minutes

24 hours working under
constant temperature

AuSn hard solder assembly
technology, higher reliability



Treatment Areas



Features

- Combined wavelength (755nm, 808nm, 1064nm) diode laser machine have a wider applicable scope.
- Fast repetition rates up to 10Hz(10 pulses-per-second), with in-motion treatment, fast hair removal for big area treatment.
- Big spot size 13x13mm², 13x39mm², high efficiency for large area treatment
- Auto alarm protection system of water flow, water level and water temperature on menu, avoids any risk at the first time
- Professional OEM/ODM services: professional services on software & hardware design and production.
- Medical CE approved, high quality promised.

Technical parameters						
Laser type	diode laser					
Display	10.4" color touch LCD screen					
Laser wavelength	808nm, 755nm, 1064nm or combined for optional					
Spot size	13*13mm ² , 13*39mm ²					
Pulse width	5~400ms adjustable					
Energy density	1-138J/CM ² , 1-150J/CM ² , 1-250J/CM ²					
Pulse frequency	1~10Hz					
Cooling mode	Water+air+Semiconductor Cooling System					
Power Supply	AC 220/50Hz ; AC 110/60Hz					
Packing	Aluminum cases					
Bar specification						
Optical Parameter						
Number of Bars	1-32					
Minimum Pitch Between Bars	1.3/2 mm					
FWHM	12x38 mm					
Electrical Parameter						
Operation Mode	CW/QCW	CW/QCW	QCW	QCW	QCW	QCW
Slope Efficiency	1.1			W/A		
Operation Temperature	25			C		
Machine specification						
Model No.	US450		US417	US415		
Output power	2700W (13x13mm ²), 3500w(13x39mm ²)		2200W(13x13mm ²)	2200W (13x13mm ²)		
Machine size	110x53x43cm		45x60x40cm	65x60x50cm		
Net weight	65.2kgs		48kgs	48.1kgs		

Micro Channel Diode Laser Machine



US460

Function


- Hair Removal
- Any hair color removal
- Any skin color hair removal
- Permanent & Painfree




US420

Micro channel diode laser :

Maximum frequency of 10HZ 808nm uniform laser diode beam penetrate through epidermis into dermis and absorbed by melanin in hair and hair follicle, producing photo-thermal effect, light can be absorbed and converted into heat damaging hair follicle tissue, so that the hair papilla is damaged without injury of surrounding tissue. Thus reaching effect of hair removal. No pain, the most safe, easy operation technology for permanent hair removal now.

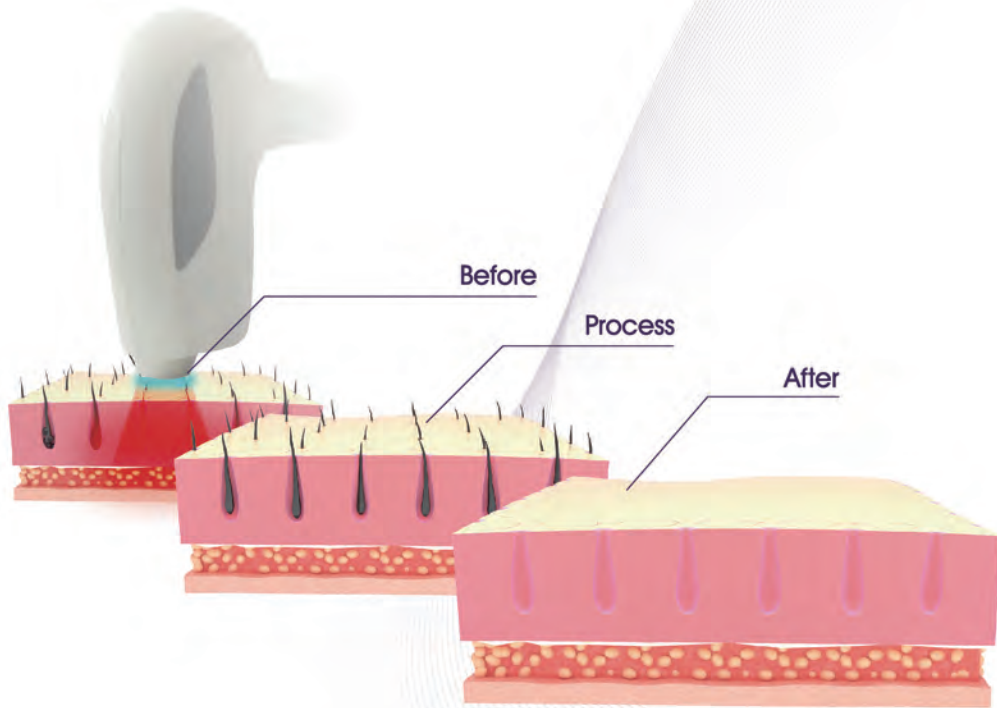

Micro channel bars increase water circulation channel inside the bars, the cooling capacity of more than 2 times the original, while the energy emitted can largely reduce the bars of the loss of their own, can achieve greater energy under the same working conditions output, to enhance the therapeutic effect.

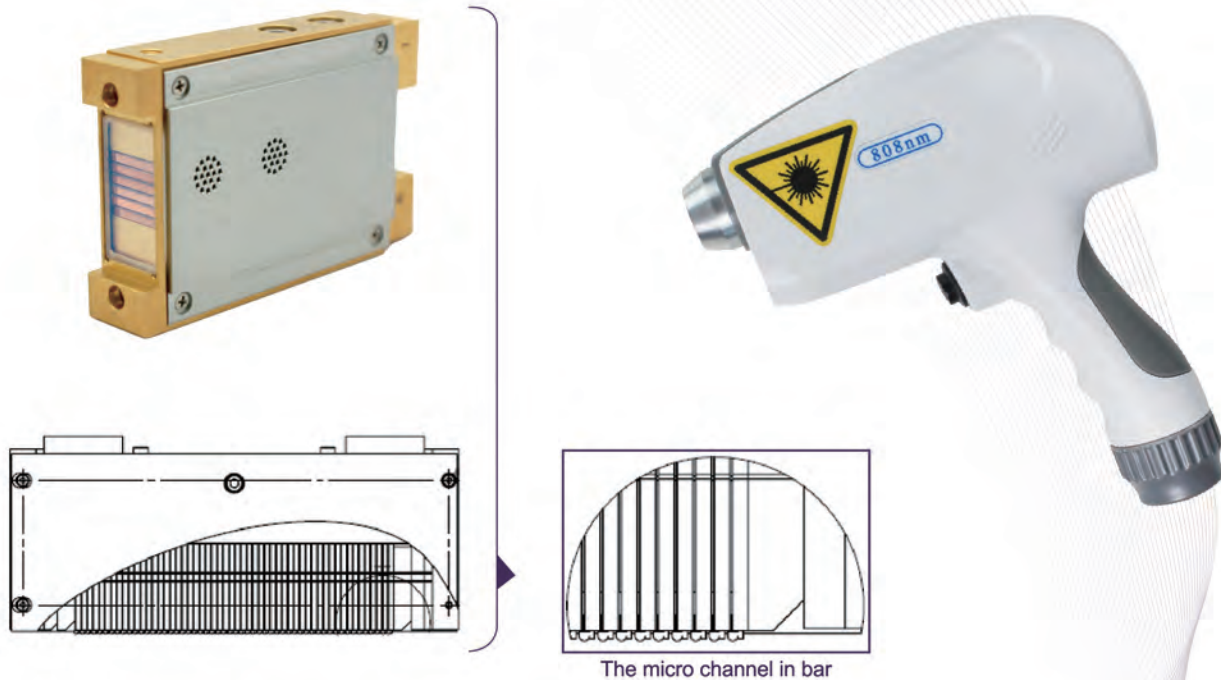

Fiber-Coupled Diode Laser

Functions

- Permanent & Painfree Hair removal on all body of skin type I-VI
- Lip hair removal, beard hair removal, chest hair removal, armpit hair removal & hair removal on outside bikini line, etc.



US417G




US408G

Handpiece



Principle

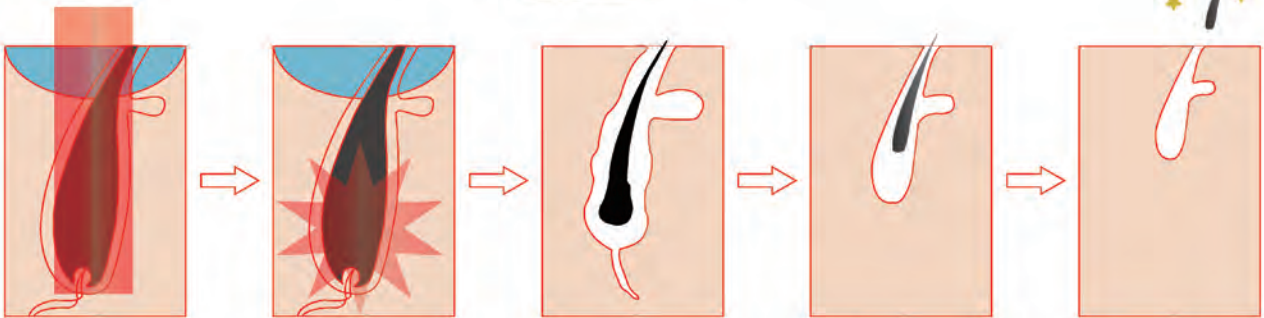

Fiber Diode Laser working principle is based on the photothermic theory. There is plenty of melanin in hair filicles and hair shafts. Melanin scatters among hair bulbs and hair shaft structure (like medullary substance, cortex and cuticula pill). Fiber Diode Laser targets on melanin and treat selectively in a precise way. Melanin could absorb the energy from the laser, which would result in temperature rapid increase, to destroy surrounding hair follicles, and finally removal hair.

Treatment pictures



Advantages

- No-channel laser generator
- No strict demand for water, no need water filter, no risk of burning laser bars.
- Biggest energy transfer
- Light handle, only 300g

Technical parameters

Technical parameters		
Model No.	US417G	US408G
Display	8" color touch screen	10" color touch screen
Laser Type	Fiber Diode Laser	Fiber Diode Laser
Wavelength	808nm	808nm
Laser Power	400W (600W for option)	600W (700W / 800W for option)
Energy	1-100J/CM ²	1-100J/CM ²
Pulse Width	5-400ms adjustable	5-500ms adjustable
Frequency	1-10HZ	1-10HZ
Spot Size	13*13mm ²	13*13mm ²
Cooling System	water+air+sapphire contact cooling	water+air+sapphire contact cooling
Crystal Temperature	0-4℃	0-4℃
Laser Water Temperature	40℃	40℃
Machine Size	54*34*38cm	99*49*46cm
Net Weight	22kg	32.2kg

IPL / SHR + Diode Laser Machine



US419

- 808nm diode laser: permanent and no pain hair removal on lip, armpit, whisker, hair line, back, chest, bikini and leg
- IPL: hair removal, skin rejuvenation, vascular therapy, pigment treatment and breast liftup.

Features

■ 808nm diode laser system:

1. Germany imported laser bars, high quality and long lifespan;
2. SHR-IN MOTION technique, fast, safe and pain free;
3. Large spot size: 13x13mm², effective results achieved;

■ 2 in 1 multifunctional beauty equipment as following:

S1: 808nm diode Laser;

S2: IPL multifunction; or SHR IPL for optional.

■ Powerful cooling system:

1. imported 316L stainless steel, better cooling, longer lifespan;

2. Semiconductor thermal management system for optional for long time working continuously.

■ IPL multifunction system:

1. 5 filters attached as 480nm, 530nm, 590nm, 640nm, 690nm;
2. 3 spot sizes for optional as 15x50mm², 12x30mm², 12x40mm²;
3. SHR IPL system for optional, 5 shots per 1s promised.

■ Germany imported ultra-quiet water pump:

1. No noise, provide a quiet and comfortable environment for patients;
2. Larger water flow, provides better cooling, and prolongs machine life.

Technical parameters						
	Controller	10.4 inch color touch LCD screen				
	Output power	3000W				
	Cooling	Water+air+semiconductor thermal management system+cooling gel				
	Handpiece No.	2 pcs standard: H1: 808nm diode laser; H2: IPL with 5 filters				
	Language	English standard, or your language to OEM				
	Package size	62x61x129cm for machine; 56x50x16 for spare parts				
	Gross Weight	78 kgs				
808nm diode laser for permanent hair removal						
	Wavelength	808nm, "Golden Standard For Permanent Hair Removal"				
	Energy Density	1-138J/cm ²				
	Therapy Mode	SHR (IN MOTION); LHR				
	Frequency	SHR: 1-10Hz; LHR: 1-3Hz				
	Pulse Width	5~400ms				
	Spot Size	13x13mm ²				
IPL for skin care therapy						
	Energy density	1-50J/cm ²		SHR IPL for skin care therapy		
	Pulse No	1~6, adjustable		1-50J/cm ²		
	Frequency	1~4s, adjustable		1~6, adjustable		
	Spot Size	12x30mm ² standard, 12x40mm ² , 15x50mm ² for optional		0.2S, 0.3S, 0.5S		
	Filters	5pcs standard: 480nm, 530nm, 590nm, 640nm, 690nm		15x50mm ² , 12x30mm ² , 8x15mm ² , for optional		
Filters						
Bars Specifications						
Optical Parameter						
	Number of Bars	1-32			1-32	
	Minimum Pitch Between Bars	1.3/2			1.3/2	
	FWHM	12x38			12x38	
	Electrical Parameters					
Operation Mode	CW/QCW	CW/QCW	QCW	QCW	QCW	QCW
Slope Efficiency	1.1			W/A		
Operation Temperature	25°C					

4H Intelligent Beauty Machine

IPL+E light+RF+YAG Laser



US002H



US609H

- E Light / IPL: hair removal, skin rejuvenation, pigment therapy, vascular therapy, breast liftup, etc.
- RF: wrinkle removal, skin tightening, body lifting, etc.
- Nd Yag Laser: tattoo removal, pigment reduction, skin rejuvenation, etc.

IPL / E light / ND YAG / RF



H1 IPL



H2 E light



H3 YAG Laser



H4 RF

Feature

- Upgradeable system, 5 options.

001: IPL System, 2 Handpieces, 5 filters;

002: E Light System, 2 Handpieces, 5 filters;

003: 3H System, 3 Handpieces as IPL, E Light, & RF;

004: 4H System, 4 Handpieces as IPL, E Light, RF, & Nd Yag;

005: SHR IPL System

- Unique designed controller, up or down- adjustable, facilitate operation, as well as save shipping cost.

- New developed luxury handpiece with streamline design & imported parts, grace, strong & durable.

- The state-of-the-art IPL filters ensures safe and effective multifunction.

- Imported water pump, no noise and larger water flow for better cooling, prolong machine life

Technical parameters

IPL system		QS Nd: Yag Laser system	
IPL energy density	1-50J/cm ²	Laser energy	1600mj
IPL frequency	1~4s	Probe	1064nm, 532nm standard (skin rejuvenation probe for optional)
Pulse no.	1~6	Instruction light	Imported infrared ray indicator
Wavelength	640nm standard	Pulse width	6ns
Spot size	15x50mm ²	Frequency	1~6Hz
E Light system		Spot diameter	4~8mm
IPL energy density	1-50J/cm ²	Model No. / US002H	
RF energy density	1-50J/cm ²	Display	8.4 inch color-touch LCD screen
IPL frequency	1~4s	Operation language	English Polish Russian Turkish German Swiss 6 languages in total
Wavelength	480nm, 530nm, 590nm, 640nm, 690nm standard	Cooling system	Japanese semiconductor +air +water + cooling gel
Spot size	12x30mm ² , 12x40mm ² optional	Output power	2500W
Bipolar RF system		Voltage	220V (110V)/15A 50Hz (60Hz)
RF energy density	1-50J/cm ²	Machine Size	44x53x89cm
RF frequency	6MHz	Package Size	57x68x107cm; 60 x30x50cm
RF delay	100-1500ms	Gross Weight	68kgs
Interval	1-4s		
Spot size	20/28/35mm dia.		

3H Intelligent Beauty Machine

IPL+E light+RF




US606S

- E Light / IPL: hair removal, skin rejuvenation, pigment therapy, vascular therapy, breast liftup, etc.
- RF: wrinkle removal, skin tightening, body lifting, etc.



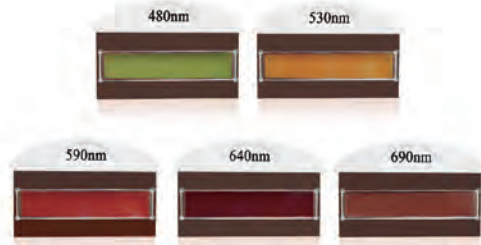

Filters

- The state-of-the-art IPL filters ensures safe and effective multifunction.
- The elegant cap facilitates inert, as well as wavelength printed to recognize for accurate treatment.
- The metal frame is produced by molding, hence, conformance and harmony.
- The filter film is pure imported; hence, the filtered wave is perfect conformity, as well as firmness.



Features

- New developed luxury handpiece with streamline design & imported parts, grace, strong & durable.
- With alarm protection system of the water flow and water temperature, the 3H beauty equipment will stop all working immediately in case of any risk at the first time, can offer the safe operation.
- USA import water connector (USA) & Germany import electronic connector
- Humanized bending design, maximum protect inner cable avoid from vertical fold, and water flow even smooth for better cooling.

Technical parameters		
IPL system		
IPL energy density	1-50J/cm ²	
IPL frequency	1~4s	
Pulse no.	1~6	
Wavelength	640nm standard	
Spot size	15x50mm ²	
E Light system		
IPL energy density	1-50J/cm ²	
RF energy density	1-50J/cm ²	
IPL frequency	1~4s	
Wavelength	480nm, 530nm, 590nm, 640nm, 690nm standard	
Spot size	12x30mm ² , 12x40mm ² optional	
Bipolar RF system		
RF energy density	1-50J/cm ²	
RF frequency	6MHz	
RF delay	100-1500ms	
Interval	1-4s	
Spot size	20/28/35mm dia.	
Model No. / US606S		
Display	7.4 inch color-touch LCD screen	
Operation language	English Polish Russian Turkish German Swiss 6 languages in total	
Cooling system	Japanese semiconductor + air + water + cooling gel	
Output power	1800W	
Voltage	220V (110V)/15A 50Hz (60Hz)	
Machine Size	43x39x36cm	
Package Size	53x54x53cm; 55x48x12cm	
Gross Weight	38kgs	

IPL/E light Beauty Machine



US002

Functions

- Hair Removal
- Skin Rejuvenation
- Pigmentation Therapy
- Vascular Therapy
- Breast Liftup, etc.



US609B

Newly designed integrated handpiece:

- Most advanced plug and play technology, no leakage of water and electricity, safe and reliable.
- USA imported water connector and Germany import electronic connector.
- Humanized bending design, 100% import tube, maximum protect inner cable from vertical fold, and water flow smoothly for better cooling.



Features

- Imported water pump, no noise and larger water flow for better cooling, prolong machine life.
- Cooling plate imported from Japan and imported sapphire ensures the best treatment effect.
- User friendly menu, maximum facilitate operation
- USA imported metal tube ensures high output power and prolongs the lifespan.
- Intelligent water flow and temperature alarm system protect machine from damage, extend machine's lifetime.

Available in 7 languages



Technical parameters	
wavelength	480nm/530nm/590nm/640nm/690nm standard, (410nm/490nm for optional)
IPL Energy Density	1-50J/cm ²
RF energy density	1-50J/cm ²
Language available	English Italian Polish Russian Turkish Spanish German 7 languages in total
Frequency	1, 2, 3, 4 s
No. of Pulse	1~6
Cooling	Japanese semiconductor + water +wind cooling +cooling gel
Voltage	220V (110V)/15A 50Hz (60Hz)
Machine specifications	
Model No	US002H US609B
Display	8.4 inch color touch screen 7.4 inch color touch screen
Output Power	2200w 1200w
Size of Spot	15x50mm ² , 12x30mm ² as standard 12x30mm ² , 12x40mm ² , 12x40mm ² , 8x40mm ² for optional 15x50mm ² for optional
Machine Size	44x53x89cm 45x43x37cm
Gross weight	60kg 34kg

Cryotherapy Slimming Machine



- Body shaping and contouring
- Body slimming and weight reduction
- Face thinner and skin tightening
- To promote and accelerate the body's metabolism
- Redundant cellulites exploded and melted

US09S



US601N

Cryotherapy

Cryotherapy is a non-invasive (without contact) procedure in which doctors place a small suction cup on a patient's skin. The device gradually extracts body heat until the subcutaneous fat is frozen.



Diode laser

The lipolysis diode laser open and release glycerol, water and free fatty acids, thus, causing to lose weight and inches.

Cavitation

Ultrasonic vibration may cause fierce impact of fat cells to make them be exploded instantaneously, reduce the amount of fat cells and thereby achieve the effects of removing fat



RF

RF can deep direct fat body, has the directional positioning RF output excellence, make fat cells organizations in fast active state, produce the cell body heat of the friction heat, in order to achieve the effect of dissolved fat.

Functions

■ Fat freeze+Vacuum+Cavitation+RF+diode laser: The treatment result will be more visible, while combine with fat freeze treatment and ultrasonic cavitation.

■ 2 fat freeze handles can work together or separately.

■ 10.4" LCD Touch Screen, beautiful and easy operation.

■ Japanese top-level Cooling plate, ensure steady cooling, good effect, long life.

Technical parameters

Output Power	1000W
Display	10.4inch color touch screen
Language	English
Dimension	59x35x112.5cm
N.W	40kg
Package	Aluminium box
Handpiece	5 pcs Handpiece / H1: Big fat freeze handle 15x5x6.5cm / H2: Small fat freeze handle 12x3.5x5cm / H3:Diode laser Big:12x9x2.6cm, Small:9.5x5.6x2.4cm / H4: Multipolar RF / H5: Cavitation
Pressure	400KPA
Temperature	-10~5degree / 30~45degree
Cooling system	Wind+water+semiconductor
Working Time	No downtime
Cavitation	Frequency: 40 KHz
Multipolar RF	Frequency:6 MHz

Cryotherapy Slimming Machine



US320C

- Body shaping and contouring
- Body slimming and weight reduction
- Face thinner and skin tightening
- Promote and accelerate body metabolism
- Redundant cellulites exploded and melted


Handpieces



Advantages

- 8 in 1 machine with many function
- 4 cryo handles can work at same time, largely save time, unique 360 degree for double chin treatment handpiece, lead the latest technology in the world
- Most comfortable treatment process: silicone rubber cryo cup with suction and release function;
- No consumable, low cost: cooling plate sealed inside cryo cup, no need anti-freezing membranes;
- Temperature from 45°C to -10°C: the automatic temperature control system provides a Heating + Cooling experience(*First Heating then cooling), which achieves greater fat removal result.
- South Korea vacuum pump imported, automatic adjustment of vacuum pressure, and ensure a quiet operation.(note: it's of big noise to the traditional vacuum pump)
- Japan double cooling plate, cooling is very stronger
- Germany casing pipe, not easily broken, strong

How it works

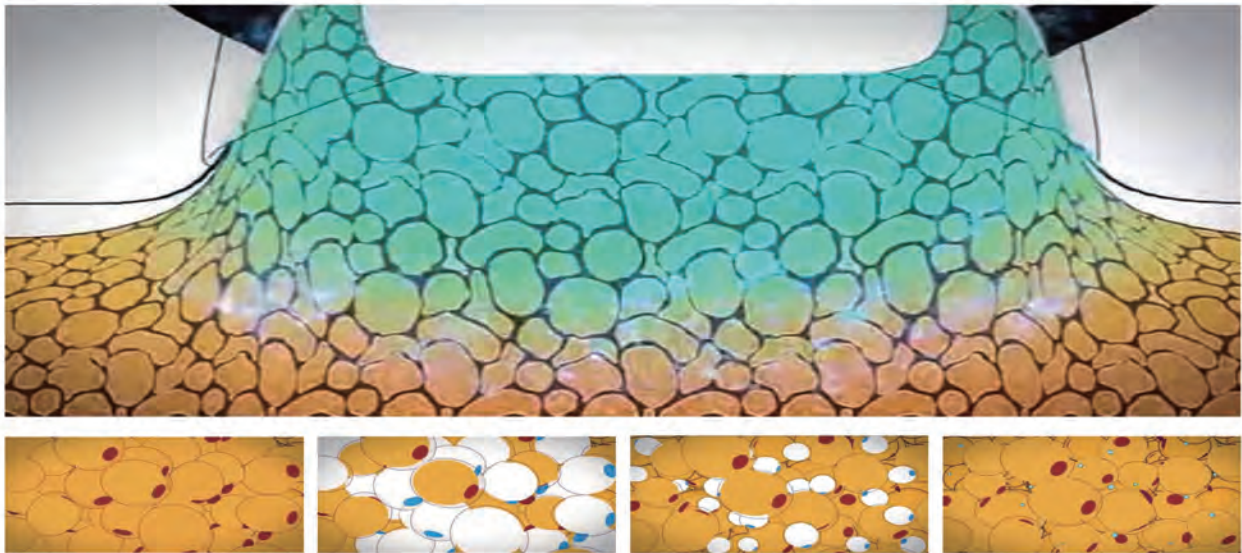

Triglyceride in fat solidifies at about 4 degrees at low temperature, cells death, decomposed through the liver, metabolism excreted out of body with urine.

Before & After



7H Beauty Slimming RF Machine



US07B



US308B

- Wrinkle removal around eyes and skin lifting;
- Body shaping;
- Neck wrinkle removal;
- Redundant cellulites exploded;
- Redundant cellulites melted;

- Body slimming and weight reduction;
- Face slimming, wrinkle removal and skin tightening;
- Fatigue relieving;
- Channels and collaterals dredging;
- Metabolism improvement and acceleration.

7H Slimming & Skin Lifting



Vacuum + RF



Multipolar RF



Tripolar RF



Bipolar RF



Cavitation



Vacuum + Roller



Microcurrent



Diode Laser

- Vacuum + RF — loss weight, skin rejuvenation, lift skin...
- Multipolar RF — Burn fat, loss weight, skin rejuvenation...
- Tripolar RF — skin rejuvenation, lift skin, care skin...
- Bipolar RF — wrinkle removal, skin rejuvenation, promote collagen regeneration...

- Cavitation — Burn fat, loss weight...
- Vacuum + Roller — loss weight, massage, metabolism
- Microcurrent — metabolism, Meridians, relieve fatigue...
- Diode laser — open and release glycerol, water and free fatty acids...

Advantages

- Integrating with 7H Beauty System, unique design in the world today;
- Vacuum pumps made in Japan;
- 8.4 inch color touch-screen, professional and easy operation;

- Unique external design, with high market competitiveness;
- No pain, no side effect.
- Multifunction, low cost and high profits;

Technical parameters

Power Supply Output	1200W	
Ultrasonic Frequency	40KHz	
Ultrasonic Power	1-50W/CM ²	
RF Frequency	6MHz	
Vacuum Pressure	400Kpa	
Roller No.	5pcs	
Current Intensity	Max 5 MA	
Pins No.	30pcs	
Shots No.	No limited	
Cooling System	Air	
Working Mode	Continuously, interval, or pulse	
Working Time	Continuous or time setting	
Model No.	US07B	US308B
Display Screen	8.4" color touch LCD screen	7.4 inch color touch screen
Language	English	English
Machine Size	40x36x112cm	58x40x33cm
Package Size	47x57x131cm	64x51x52cm
Gross Weight	46Kg	23kg
Voltage	220V/110V	220V/110V

Vacuum Slimming Machine

Application

Face:

Smooth wrinkles, firm the skin, clarify the complexion, reduce double chin, black circle around eyes removal, upper eyelid lifting and wrinkle removal, all face treatment.

Body:

Smooth cellulite, firm the skin, release fat, lighten legs, body contouring after birth, stretch mark improvement, body shaping and slimming , all body treatment.



US08S

Handpieces




- **Vacuum negative pressure** can suck the skin to fill the blood vessels. This can speed up blood circulation and promote metabolism. In addition, vacuum negative pressure can bring a lot of lipolytic enzymes, these decomposition enzymes in the role of radio frequency can be better to dissolve fat.
- **Mechanical roller**
Mechanical roller can promote lymphatic reflux and accelerate the metabolism of mind, so that the resulting physiological waste is quickly excreted.



- **20Mhz RF**, three sizes treatment heads.
 - Small size head** — Utilized RF & integrated skin cooling mechanism to treat lax, sagging skin on the face and body.
 - Middle size head** — Utilized RF and vacuum, for cellulite reduction of localized fatty deposits in small and sensitive areas such as the chin, neck, inner thighs and under arm.
 - Big size head** — Utilized RF and vacuum for cellulite reduction in large areas such as the thighs, buttocks, arms and abdomen.

Principle

Acting simultaneously on the release of stubborn localized fat and on skin quality (firmness, cellulite aspect), the new tech (Roller and motorized flap combined with sequential suction) allows the body technology to act directly on the adipocytes present in our hypodermis (slimming cells) and on the fibroblasts present in our dermis (rejuvenating cells)

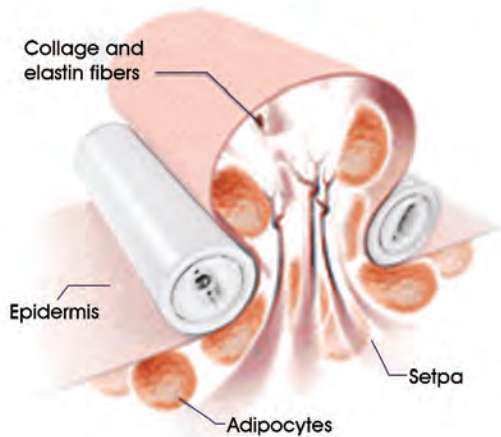

3 in 1 CO2 Fractional Laser Machine

Fractional, Normal, V²Lr



Functions

- Skin renewing and resurfacing
- Wrinkle removal, skin tightening
- Acne and acne scar removal
- Smooth burnt scars and surgery scars, etc
- Remove intractable chloasmas and pigmentation
- Sun damage recovery.
- Vaginal tightening laser:
 1. Contracting the vagina
 2. Labium Whitening
 3. Solve the problem of dyspareunia and Colpoxerosis
 4. Young vaginal tissue

US800



US900


Work principle

Co2 fractional laser: use 10.6 μ m invisible light, by the treatment method of non-contact cut, burn, and exposure, the heat will congeal, ablative vaporized, or carbonize the tissues, so as to shorten the operation time, and patients will have less bleeding, healing fast, and greatly enhance anti infectious ability.



Vaginal tightening

It concentrated thermal heating of the vaginal tissue inner layers induces collagen and elastin contraction and regeneration of these proteins in the longer term. These molecular signals activate the entire wound healing process of the body's local immune system.

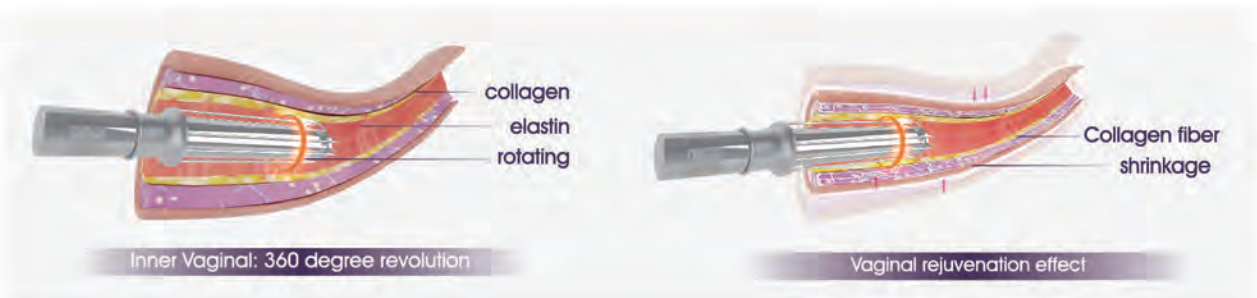

Fractional: 7 Scan Graphics



Normal Co2 : 3 Working Modes

		
Single	CW	Pulse

3 in 1 system, easy and convenient for operation.



fractional CO2 laser



normal CO2 laser



vaginal probe 1



vaginal probe 2



vaginal probe 3

Features

■ USA Coherent RF laser emitter

1. Even output energy, deeper penetration (8-10mm under the skin)
2. Long lifetime: 8-12 years
3. Smaller spot size: 0.01mm at dia

■ 6 scan modes: sequence, random, maximum distance A&B, Equidistance and fast shot;

■ Imported excellent 7 articular optical-arm, easily operated and greatly reduces the energy loss;

■ Ultra-quiet pump used, lower noise but larger water flow;

■ 3 in 1 system, easy and convenient for operation.

1. Fractional CO2 Laser for skin rejuvenation and scar removal
2. Normal CO2 laser for surgery operation
3. Vaginal tightening for optional.


■ 7 scan shapes: square, rectangle, round, triangle, oval, 6-diamond shape and line

■ The screen can be folded up and down, rotate 180° left and right

Technical parameters		
Wavelength	10.6μm, far-infrared laser	
Pulse frequency	33.3 Hz	
Duration	0.1~10ms	
Distance	0.1~2.6mm	
Interval	0~6s	
Repeat	1~5000ms	
Overlap	1~20 times	
Working Mode	Continuous pulse, single pulse and interval pulse	
Scan Graphics	Sequence scan, Random scan, Maximum distance A&B, Equidistance, and Fast shot	
Pulse energy	10mj to 130mj is optional for each dot.	
Display	10.4 inch color touch screen	
Machine specifications		
Model No	US800	US900
Laser emitter	USA Coherent RF laser emitter	glass tube laser emitter
Spot size	0.01mm at dia	0.2mm at dia
Average power	50W	30w
Cooling system	Air	water+air
Size of machine	58.5cmx43.8cmx107cm	72cmx38.5cmx42.5cm

3 in 1 CO2 Fractional Laser Machine


Fractional, Normal, V²Lr

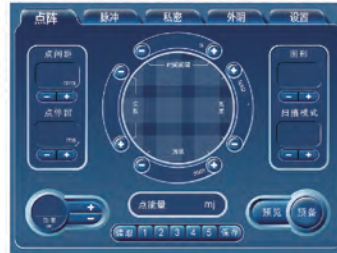

- 3 in 1 system:
 1. Fractional CO2 laser;
 2. Normal laser for surgery ;
 3. Individual Gynae and Vulva .



- Freely draw graphics as required, providing drawing pen.

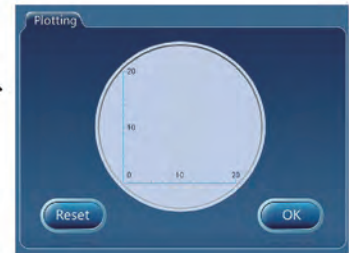

US800N



US900N



- 3 tips for CO2 fractional laser, with different penetrating depth.

Picosecond Laser Machine



US506

- All color tattoo, eyeliner and lip line removal
- Epidermis and dermis pigment
- Nevus of Ota, blue nevus, black nevus, coffee spot

- Senile plaques, freckle, sunburn spot and chloasma
- Skin whitening and wrinkle removal etc

Treatment Tips


US506A



US506B



Theory

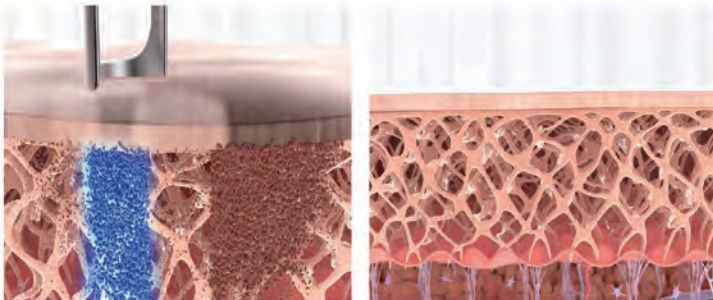

- Picosecond laser machine uses a unique wavelength of light and a specialized lens to convert laser energy into gentle pressure. This pressure squeezes cells, activating the natural cell signaling processes that create new collagen and elastin. Picosecond laser will minimize the side effects of thermal effect, without burning or damaging the skin, can achieve the goal of almost solve all kinds of pigment spots, is better than traditional laser spot whitening effect.

Advantages

- 1. High stability(our patent and advantage): with advanced full-automatic water temperature protection cooling system;
- 2. Elegant metal shell, strong and durable;
- 3. Fast Repetition: up to 15HZ, quick and efficient treatment;
- 4. High output energy: up to 800MJ, most effectiveness;
- 5. Tips including carbon peel and honeycomb tip.



Technical parameters

Laser Type	Picosecond laser	
Spot diameter	2-10mm	
Cooling type	Water circulation inside + air cooling	
Model	US506A	US506B
Frequency	1-10Hz	1-15Hz
Pulse width	≤6ns	≤6ns
Pulse energy	1064nm-800mJ 532 nm-500mJ 755nm-500mJ	1064nm-800mJ 532 nm-500mJ 755nm-500mJ
POWER	600W	800W
Machine size	36*55*50cm	36*55*50cm
Dimensions	49*59*67cm	49*59*67cm
G.W	36.7KG	36.7KG

Picosecond Laser Machine



US508

- All color tattoo, eyeliner and lip line removal
- Epidermis and dermis pigment
- Nevus of Ota, blue nevus, black nevus, coffee spot

- Senile plaques, freckle, sunburn spot and chloasma
- Skin whitening and wrinkle removal etc

Treatment Tips




Normal Tip



Honeycomb Tip

Advantage

- PicoSecond laser machine used unique Honeycomb Focused technology to form a skin effect vacuolization, which can protect skin from damage during treatment.
- PicoSecond laser machine make tattoo & pigment removal treatment process from 5 to 10 times reduced to 2-4 times, greatly reduce the treatment and recovery time, with fast and obvious effective.
- It can removal all kinds of pigment and tattoo effectively and safely, because picosecond laser use precise positioning of the target tissue therapy to reduce damage to the skin to achieve the freckle effect.
- The picosecond laser use ultra-short pulses (one trillionth of a second in length) to hit the melanin with a great pressure, the melanin shatters into tiny dust-like particles. Because the particles are so small, they are more readily absorbed and eliminated by the body. It will significantly reduce postoperative swelling, melanin precipitate phenomenon.

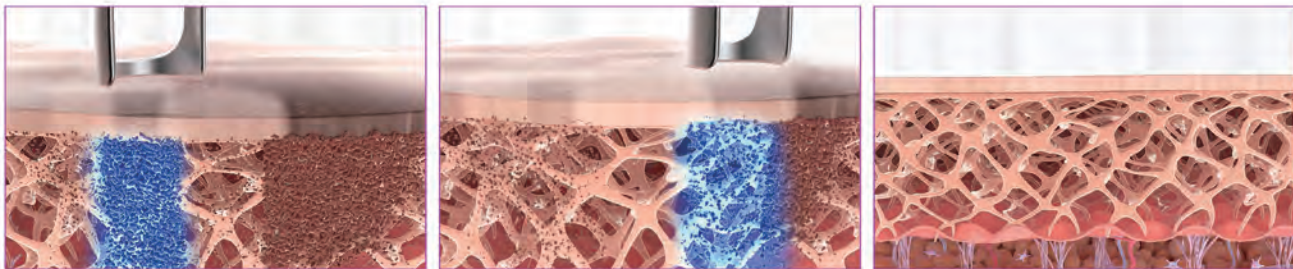

Start treatment

During treatment

Finish the treatment

Technical parameters

Technical parameters			
Laser Type	Picosecond laser		
Spot diameter	2-10mm		
Cooling type	Water circulation inside + air cooling		
Mode	US508A	US508B	US508C
Frequency	1-10Hz	1-10Hz	1-10Hz
Pulse width	800PS	500ps	300ps
Pulse energy	1064nm-1200mJ 532nm-600mJ	1064nm-800mJ 532 nm-400mJ	1064nm-500mJ 532 nm-250mJ
POWER	1200W	2000W	2000W
Machine size	50*85*160cm	54*88*150cm	54*88*150cm

Q-Switch ND: YAG Laser Machine



US407



US412

- Pigment deposit dispelling
- Skin Rejuvenation
- Tattoo removal: Remove black, blue, brown, and colorful (such as red, yellow, green) tattoo on eyebrow, eye line, lip line, and other parts of body
- Eliminates the age pigments, spots, birthmarks, and pigment changes

Treatment head & Electro optic Q-switched

Electro-optic Q-switch ND yag laser, releases actual single pulse of laser energy.




Normal Probe



Skin Rejuvenation Probe



Q-switched

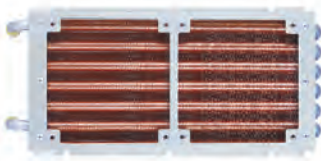

Stainless steel Radiator

Imported 316L stainless steel, better cooling, longer lifespan

Imported water pump

Germany imported water pump, no noise and larger water flow for better cooling, prolong laser life and enjoy treatment in salon.



Features

- Electro-optic Q-Switch ND:YAG Laser, achieve virtually painless treatment and long lasting result
- 3 Working Modes provides more extensive treatment applications
- Advanced adjustable Aiming Beam assure accurate treatment even on tiny pores.
- Spot size: 2~10mm adjustable
- Alarm protection system of water flow and water temperature: protect people and machine against any risk at the first time.
- Korea imported Laser Arm: stable laser output, and maximum reduces the loss of laser energy

Technical parameters

Laser type	Q-Switch ND : YAG Laser	
Wavelength	1064nm/532nm	
Spot diameter	2-10 mm, with adjustable treatment probes	
Laser aiming	Red semiconductor aiming light, wavelength is 650-670nm	
Pulse width	2-6 ns @ Q-Switch Mode; 300us @ Quasi-Long Pulse	
Frequency	1-10Hz	
Cooling system	Closed-off water circulation inside + air cooling	
Voltage	220V/110V	
Model No.	US407	US412
Power	1200W	1500W
Single pulse energy	1064nm-600mJ; 532nm-300mJ	1064nm-800mJ; 532nm-400mJ
Double pulse energy	1064nm-1200mJ; 532nm-400mJ	1064nm-1600mJ; 532nm-500mJ
Quasi-long pulse energy	1064nm-2J	1064nm-3J
Machine size	44*70*170 cm	38*88*185cm
Gross weight	104Kgs	115kg

ND YAG Laser Machine



Functions

- Pigment deposit dispelling
- Skin Rejuvenation
- Tattoo removal: Remove black, blue, brown, and colorful (such as red, yellow, green) tattoo on eyebrow, eye line, lip line, and other parts of body
- Eliminates the age pigments, spots, birthmarks, and pigment changes

US409



US406

Treatment heads

Two wavelengths: 1064 nm, 532 nm



Treatment tips for US409



Treatment tips for US406

Treatment pictures



Features

- Fashionable design
- Big laser output energy: higher efficiency and more comfortable treatment
- Humanized treatment menu: English system language, simple operation
- Alarm of water flow: if no water inside or few water inside machine, the system will alarm automatically-around alarm and stop all working immediately

- 100% American imported "plug and play" handle connector, integrating with a perfect water-electricity isolation assembly inside; increase the stability of the machine greatly and facility actual operation actually
- High quality machine shell of ABS material, with OEM painting services

Technical parameters

Wavelength	1064/532nm	
Screen	7.4" color-touch screen	
Joint part	Adopts the most advanced (Plug-and-play) joint part	
Type of laser	sapphire and ruby switch Q/KTP/YAG laser instrument	
Pulse energy	1600mJ	
Instruction light	infrared ray indicator	
Width of pulse	6ns	
Frequency	1 to 6Hz	
Spot diameter	1-8mm	
Cooling system	Wind + water	
Voltage	220V(110V)/5A 50Hz	
Machine specifications		
Model No.	US406	US409
Size	55x40x39cm	38x52x110cm
Weight:	15kg	56kg



ND YAG Laser Machine



- 4 different tips
(1) 1064nm
(2) 532nm
(3) Carbon peeling
(4) Skin rejuvenation

- User Friendly Platform, Easy Operation



- Water Alarm System, prolong machine lifespan

- Aiming beam brightness adjustable assure accurate treatment




US400N

Diode Laser Vascular Removal Machine

Functions

- Vascular lesions
- Proliferative lesion
- Age spots and sun spots
- Spider veins
- Linear angiectasis
- Cherry angiomas

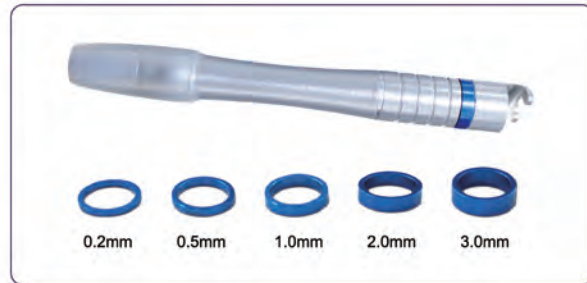

980nm diode laser with infrared aiming beam be easy to treatment



US410C

Specification

Laser type	diode laser
Wavelength	980nm +/- 10nm
Display	8.4" color touch screen
Output power	30W
Spot size	0.4-0.5mm
Pulse width	1-100ms
Frequency	1-10Hz
Time setting	1-10minutes
Aiming beam	Red semiconductor aiming light
Cooling system	Air cooling
Voltage	220V 60Hz / 110V 50Hz
Dimensions	40*38.5*20cm
Net weight	8kgs
Package size	50*40*48cm
Gross weight	17kgs

H₂O₂ Facial Jet Peel Machine



Functions

- Deep cleaning
- Eye pouch improvement
- Skin texture improvement
- Allergic skin improvement
- Water replenishing to skin
- Skin fatigue improvement
- Black rim of eye improvement
- Acne removal/acne scar improvement

US710A




US710HM


How it works

Oxygen H₂O₂


The modern medical science thinks that The decay is the process of acidification of matter (oxidation).absorb O₂, drinking and smoking, environmental pollution etc could produce oxygen free radicals in human body. This could vandalize the cell tissue, Genetic disease and aging in the body. Hydrogen could effectively remove free radicals in the body.Hydrogen anti-oxidation is much more powerful than vitamin C, carrot, lecithin etc people already familiar antioxidant substance.

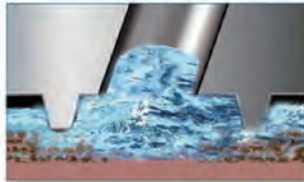

Step 1 **Cleansing and Exfoliation**
Dead skin cells are removed to reveal healthy new skin.

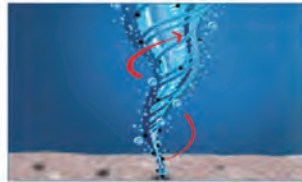

Step 2 **Acid Peel**
This gentle peel helps loosen dirt and debris from pores without irritation.

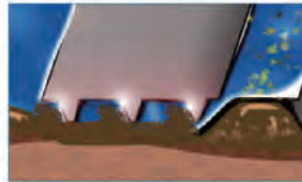

Step 3 **Extractions**
Painless, automated extractions use vacuum suction to draw out impurities.



Step 4 **Hydration**
Antioxidants and Hyaluronic Acid are Vortex-Fused™ to nourish and protect.

Handles



● Bubble Cleaner Handle



● Spray Handle



● Ultrasound Handle



● Skin Scrubber Handle



● RF Handle



● Cold Handle



7 Colors 149 LEDs



● Mask

Technical parameters

Model No.	710HM
Name	H ₂ O ₂ Facial Jet Peel Machine
Display	8" color touch screen
Vacuum Range	80Kpa
Radio Frequency	1Mhz, Bipolar
Lon Lifting	500Hz(Digital Ion Lifting)
Ultrasound	1Mhz 12W/cm ²
Machine Size	46*43*30cm
Net Weight	20kg
Voltage	220V / 110V

Air Cooling Machine

Functions

- Can work with 808nm diode laser hair removal machine,
- CO2 fractional laser machine,
- SHR machine, IPL machine and Nd yag laser machine ...
- Suitable for all machines if patients feel pain and needs cooling.
- Offer a much comfortable treatment, say goodbye to pain.



Cool-01

Technical parameters

Technology	Air cooling Machine
Power	800W
Air volume	0—280 L/Min
Voltage	200V
Lowest temperature	—35°C
The length of pipe	180cm
Machine size	43x68x77cm
Weight	55KGS



CERTIFICATES

To ensure that our products are truly world class, safe and effective, all GlobalIPL products are designed and manufactured strictly in accordance with ISO:13485 and are certified by CE . Please see below our certifications.



EPILATION
808nm Diode Laser
Dye Pulsed Light
SHR/IPL
808nm Fiber Laser



GYNECOLOGICAL
Vaginal Tightening
Vaginal Whitening
Urinary Incontinence




**ANTI-AGING AND
CELLULITE**
Radio Frequency
Co2 Fractional Laser
Vacuum+RF+Infrared



SKIN IMPERFECTION
Radio Frequency
Dye Pulsed Light
SHR/IPL
Co2 Fractional Laser
Q-Switch Laser/Pico laser

GlobalIPL factory located in Dezhou, Shandong province. It has over 7,000 square meters semi-automated production workshop with a monthly production capacity of 10,00 units, more than 200 employees. We could assure you the timely delivery ,3-5 working days.

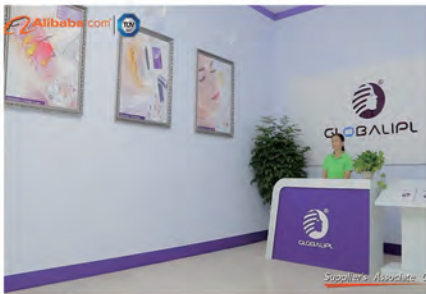

Entrance



Showroom



Production line




Warehouse



QC department



Packing

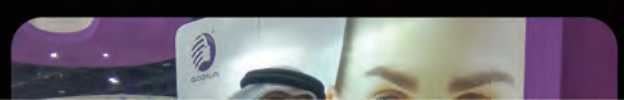

EXHIBITION

DUBAI

ITALY

SPAIN

HONG KONG

MILAN

TURKEY

.....

.....